

Confidence without compromise.

The TASKA CX™ is the hand destined to play a big role in users' daily lives.

Ground breaking motor design

Perform tasks faster, with more precision and superior grip strength.

Capacitive touch

All fingertips, including the thumb, have capacitive touch and can be used on any common touch screen.

One-click Grip Tip replacement

Grip Tips can be easily replaced with a single simple click.

Intuitive button layout

Cycle through grips, control Bluetooth™ and more with the ergonomically designed button panel.

Coverplate customisation

Replaceable Coverplates are available in a range of colours and designs.

Anatomical design

Contoured lines and discreet Wrist Flexion buttons create a streamlined profile.



For more information, visit taskaprosthetics.com

The **most advanced** prosthetic hand on the market.

The **TASKA CX** makes daily tasks easier, giving the user an enhanced sense of connection via fine control, contoured aesthetics, and personalisation features.

Robust and versatile

- Shock absorbing Knuckle Block.
- Resettable Knuckle Release for reliability and durability.
- Perform a wide range of daily activities.
- All the robustness you've come to expect from TASKA.

Waterproof

- No need for fragile gloves due to IP67 rating.
- Wash both hands easily and confidently with soap and water.
- Complete everyday tasks like washing dishes, cleaning the car, etc.
- Wrist option available with TASKA HydroSeal™ Waterproofing Technology allowing submersion past the wrist.

Easy

- Finger splay for reliable grips, lock button for safety and control, motorised thumb, and customisable grips using the MyTASKA app.
- Intuitive control systems, pattern recognition compatible, and improved button layout.
- Inbuilt Flexion Wrist.

Options

- Quick Disconnect Wrist (QDW) or Low Profile (LP).
- Available in small and medium sizes.
- Black or sand base.
- A range of Coverplates are available.

Warranty and service


- 2-year standard warranty (extendable to 5-years).
- Complimentary annual service.
- Loan units provided during all servicing and repairs.

Package contents

- Spare Coverplate & replacement Grip Tips.
- Lamination Collar Kit (Low Profile only).

All TASKA hands are compatible with a wide variety of myoelectric components. See the TASKA website for full details.

Dimensions and weights

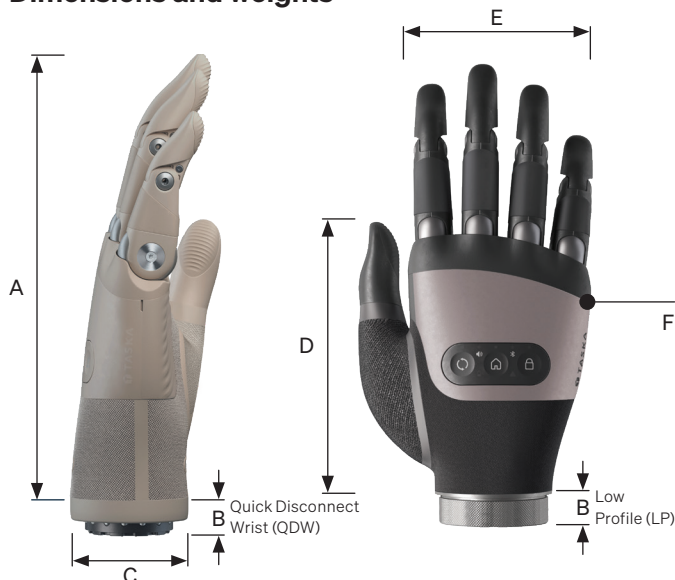

DIMENSIONS			SMALL	MEDIUM
A	Hand length	(QDW) (LP)	172 mm 168 mm	180 mm 176 mm
B	Hand base to wrist base	(QDW) (LP)	16 mm 12.5 mm	16 mm 12.5 mm
C	Diameter at wrist		48 mm	54 mm
D	Thumb tip length	(QDW) (LP)	112 mm 108 mm	114 mm 110 mm
E	Palm width		76 mm	83 mm
F	Knuckle circumference		184 mm (7¼")	197 mm (7¾")

KEY SPECIFICATIONS TABLE			
Maximum carry load		20 kg	
Grip forces	Thumb	24 N	
	Finger	13 N	
Digit speed	Thumb	0.52 s or 120°/s	
	Finger	0.7 s or 120°/s	
	Thumb rotator	0.55 s or 160°/s	
Maximum opening	Small	100 mm (3.94")	
	Medium	110 mm (4.33")	
Wrist rotation range	(QDW)	Unlimited	
	(LP)	300°	
Wrist extension lock position		15°	
Wrist Flexion lock position		25°	
Waterproof rating		IP67	
Operating temperature range		-10°C to 45°C (14°F to 113°F)	
Operating voltage		7.4 V (nominal)	
Weight	Small	(QDW)	471 g
		(LP)	430 g
	Medium	(QDW)	533 g
		(LP)	481 g