

Confidence without compromise.

TASKA CX™

Fully submersible

Using TASKA's Hydroseal Waterproofing Technology.

Touch screen compatible

Capacitive touch across all digits can be used on any common touch screen.

Damage protection

Patented shock absorbing Knuckle Block designed to handle day-to-day knocks.

Anatomical design

With natural proportions, anthropomorphic curves and smooth, quiet movements.

Intuitive control

Grip changes and locking accessible via MyoTriggers or the button panel.

Built-in wrist flexion

The world's lightest multi-articulating hand with adjustable wrist flexion and extension.



do more™

 **TASKA™**
by eqwal.

The most advanced prosthetic hand on the market.

The TASKA CX combines **refined** design with **reliable** performance. Fully waterproof and multi-articulating, it features powered thumb rotation, integrated wrist flexion, and a lightweight build for an **empowering** experience every day.

Options

- Quick Disconnect Wrist (QDW) or Low Profile (LP).
- Available in small and medium sizes.
- Black and sand colours available.
- A range of Coverplates are available.

Compatibility

- TASKA PowerCore system, Lamination Collars, Seal Rings and TASKA Gloves.
- A wide variety of myoelectric components. See taskaprosthetics.com for full details.

Warranty and service

- 2-year standard warranty (extendable to 5-years).
- Complimentary annual service.
- Loan units provided during all servicing and repairs.

Package contents


- Spare Coverplate & replacement Grip Tips.
- Lamination Collar Kit (Low Profile only).

Support

- The MyTASKA app can be used to set up, customise and train on any TASKA hand.
- Educational support is available through the TASKA Learning Academy.
- In-market support is provided for the lifetime of the product.

To learn more talk to your local TASKA representative or visit taskaprosthetics.com

Dimensions and weights

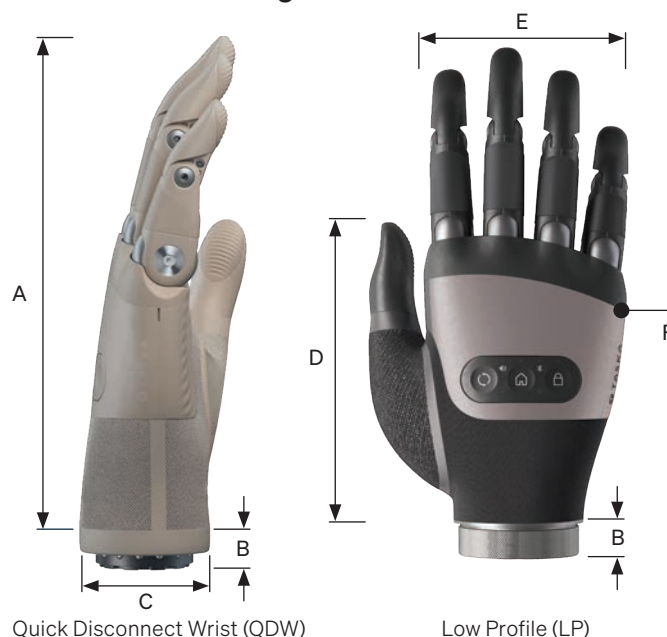

DIMENSIONS			SMALL	MEDIUM
A	Hand length	(QDW) (LP)	172 mm 168 mm	180 mm 176 mm
B	Hand base to wrist base	(QDW) (LP)	16 mm 12.5 mm	16 mm 12.5 mm
C	Diameter at wrist		48 mm	54 mm
D	Thumb tip length	(QDW) (LP)	112 mm 108 mm	114 mm 110 mm
E	Palm width		76 mm	83 mm
F	Knuckle circumference		184 mm (7¼")	197 mm (7¾")

KEY SPECIFICATIONS TABLE			
Maximum carry load			25 kg
Grip forces	Thumb		Up to 29 N
	Finger		Up to 23 N
Digit speed	Thumb		0.52 s or 120°/s
	Finger		0.7 s or 120°/s
	Thumb rotator		0.55 s or 160°/s
Maximum opening	Small		100 mm (3.94")
	Medium		110 mm (4.33")
Wrist rotation range	(QDW)		Unlimited
	(LP)		300°
Wrist extension lock position			15°
Wrist Flexion lock position			25°
Waterproof rating			IP67
Operating temperature range			-10°C to 45°C (14°F to 113°F)
Operating voltage			7.4 V (nominal)
Weight	Small	(QDW)	471 g
		(LP)	430 g
	Medium	(QDW)	533 g
		(LP)	481 g

do more™

T TASKA™
by eqwal.